

## Caspase 6 (Mch 2) Ab-3 (Clone 6CSP03; same as 14.1.190)

### Mouse Monoclonal Antibody

Cat. #MS-1126-P0, -P1, or -P (0.1ml, 0.5ml, or 1.0ml at 200µg/ml) (Purified Ab with BSA and Azide)

Cat. #MS-1126-P1ABX or -PABX (0.1ml or 0.2ml at 1.0mg/ml) (Purified Ab without BSA and Azide)

**Description:** Caspases play a major role in the transduction of the apoptotic signal and execution of apoptosis in mammalian cells. Caspase 3 (CPP32) and Caspase 6 (Mch2) are the major active caspases in apoptotic cells, and are activated in response to distinct apoptosis-inducing stimuli and in all cell lines analyzed. Both CPP32 and Mch2 are present in apoptotic cells as multiple active species. Caspase-6 cleaves nuclear mitotic apparatus protein (NuMA) and mediates the shrinkage and fragmentation of nuclei.

**Mol. Wt. of Antigen:** 34kDa

**Epitope:** Caspase 6 prodomain

**Species Reactivity:** Human. Others-not tested.

**Clone Designation:** 6CSP03 (same as 14.1.190)

**Ig Isotype / Light Chain:** IgG<sub>1</sub> / κ

**Immunogen:** Recombinant human caspase 6 prodomain.

### Applications and Suggested Dilutions:

- Western Blotting (Ab 1µg/ml for 2hrs at RT)
- Immunoprecipitation (Denatured verified)  
(Use Protein G) (Ab 2µg/mg protein lysate)
- Immunohistology (Formalin/paraffin)  
(Use Ab at 1-2µg/ml for 30 min at RT)
- \* [Staining of formalin-fixed tissues REQUIRES boiling tissue sections in 10mM citrate buffer, pH 6.0, (**NEOMARKERS'** Cat. #AP-9003), for 10-20 min followed by cooling at RT for 20 min.]

The optimal dilution for a specific application should be determined by the investigator.

**Positive Control:** Jurkat cells. Tonsil.

**Cellular Localization:** Cytoplasmic

**Storage and Stability:** Ab with sodium azide is stable for 24 months when stored at 2-8°C. Antibody WITHOUT sodium azide is stable for 36 months when stored at below 0°C.

### Supplied As:

200µg/ml antibody purified from the ascites fluid by Protein G chromatography. Prepared in 10mM PBS, pH 7.4, with 0.2% BSA and 0.09% sodium azide. Also available without BSA and azide at 1mg/ml.

### Suggested References:

1. Fernandes-Alnemri T, et. al. (1995) Cancer Res, 55:2737-2742.

### Limitations and Warranty:

Our products are intended FOR RESEARCH USE ONLY and are not approved for clinical diagnosis, drug use or therapeutic procedures. No products are to be construed as a recommendation for use in violation of any patents. We make no representations, warranties or assurances as to the accuracy or completeness of information provided on our data sheets and website. Our warranty is limited to the actual price paid for the product. NeoMarkers is not liable for any property damage, personal injury, time or effort or economic loss caused by our products.

### Material Safety Data:


This product is not licensed or approved for administration to humans or to animals other than the experimental animals. Standard Laboratory Practices should be followed when handling this material. The chemical, physical, and toxicological properties of this material have not been thoroughly investigated. Appropriate measures should be taken to avoid skin and eye contact, inhalation, and ingestion. The material contains 0.09% sodium azide as a preservative. Although the quantity of azide is very small, appropriate care should be taken when handling this material as indicated above. The National Institute of Occupational Safety and Health has issued a bulletin citing the potential explosion hazard due to the reaction of sodium azide with copper, lead, brass, or solder in the plumbing systems. Sodium azide forms hydrazoic acid in acidic conditions and should be discarded in a large volume of running water to avoid deposits forming in metal drainage pipes.

### For Research Use Only

Thermo Fisher Scientific  
Anatomical Pathology  
47777 Warm Springs Blvd.  
Fremont, CA 94539, USA  
Tel: 1-510-991-2800  
Fax: 1-510-991-2826  
<http://www.thermo.com/labvision>



Manufactured by:  
NeoMarkers  
For  
Lab Vision Corporation

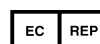

Thermo Fisher Scientific  
Anatomical Pathology  
93-96 Chadwick Road, Astmoor  
Runcorn, Cheshire WA7 1PR, UK  
Tel: 44-1928-562600  
Fax: 44-1928-562627  
[Labvision.uk@thermofisher.com](mailto:Labvision.uk@thermofisher.com)